


## Bioengineering 280A Principles of Biomedical Imaging


Fall Quarter 2007  
MRI Lecture 3

Thomas Liu, BE280A, UCSD, Fall 2007

## Fourier Sampling

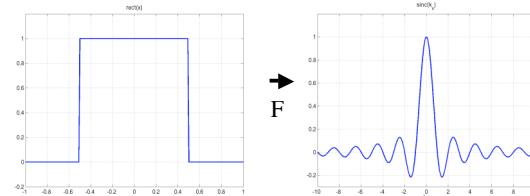


Instead of sampling the signal, we sample its Fourier Transform

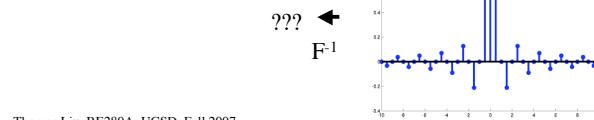

Thomas Liu, BE280A, UCSD, Fall 2007


## Fourier Sampling

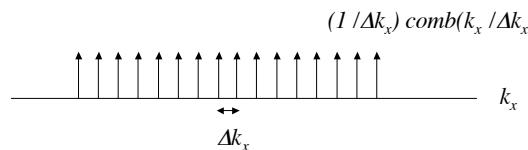

$$\begin{aligned} G_S(k_x) &= G(k_x) \frac{1}{\Delta k_x} \text{comb}\left(\frac{k_x}{\Delta k_x}\right) \\ &= G(k_x) \sum_{n=-\infty}^{\infty} \delta(k_x - n\Delta k_x) \\ &= \sum_{n=-\infty}^{\infty} G(n\Delta k_x) \delta(k_x - n\Delta k_x) \end{aligned}$$

Thomas Liu, BE280A, UCSD, Fall 2007

## Fourier Sampling -- Inverse Transform

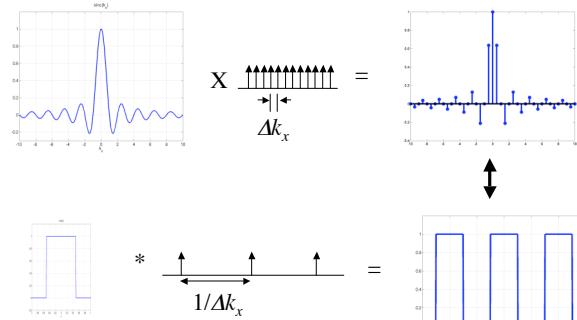


Thomas Liu, BE280A, UCSD, Fall 2007

## Fourier Sampling -- Inverse Transform

$$\begin{aligned}
 g_s(x) &= F^{-1}[G_s(k_x)] \\
 &= F^{-1}\left[G(k_x) \frac{1}{\Delta k_x} \text{comb}\left(\frac{k_x}{\Delta k_x}\right)\right] \\
 &= F^{-1}[G(k_x)] * F^{-1}\left[\frac{1}{\Delta k_x} \text{comb}\left(\frac{k_x}{\Delta k_x}\right)\right] \\
 &= g(x) * \text{comb}(x\Delta k_x) \\
 &= g(x) * \sum_{n=-\infty}^{\infty} \delta(x\Delta k_x - n) \\
 &= g(x) * \frac{1}{\Delta k_x} \sum_{n=-\infty}^{\infty} \delta(x - \frac{n}{\Delta k_x}) \\
 &= \frac{1}{\Delta k_x} \sum_{n=-\infty}^{\infty} g(x - \frac{n}{\Delta k_x})
 \end{aligned}$$

Thomas Liu, BE280A, UCSD, Fall 2007


## Nyquist Condition

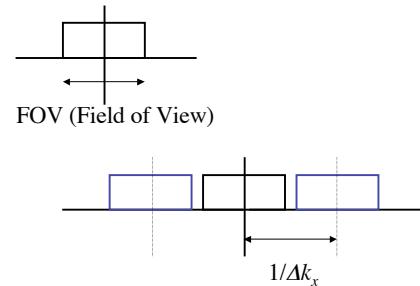

To avoid overlap,  $1/\Delta k_x > \text{FOV}$ , or equivalently,  $\Delta k_x < 1/\text{FOV}$

Thomas Liu, BE280A, UCSD, Fall 2007


## Aliasing

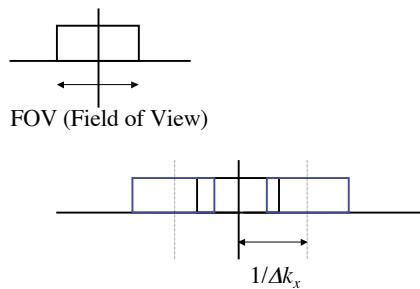





Aliasing occurs when  $1/\Delta k_x < \text{FOV}$


Thomas Liu, BE280A, UCSD, Fall 2007

## Intuitive view of Aliasing

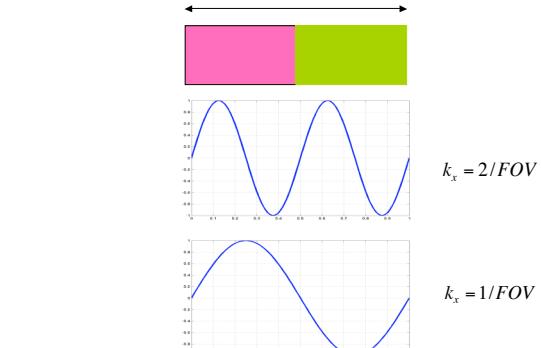

Thomas Liu, BE280A, UCSD, Fall 2007

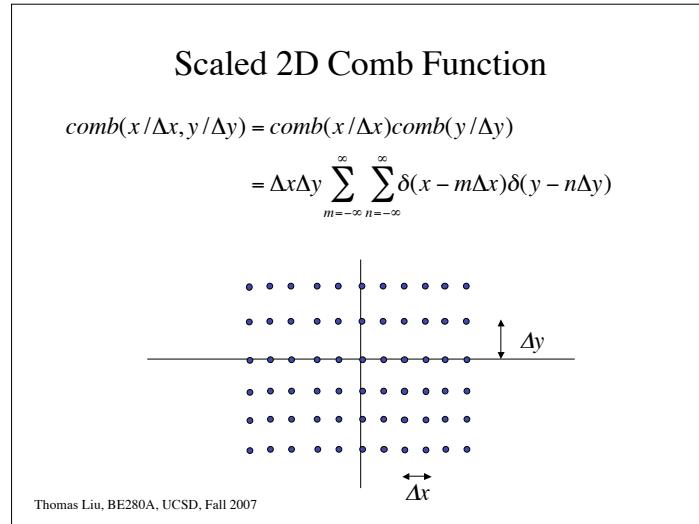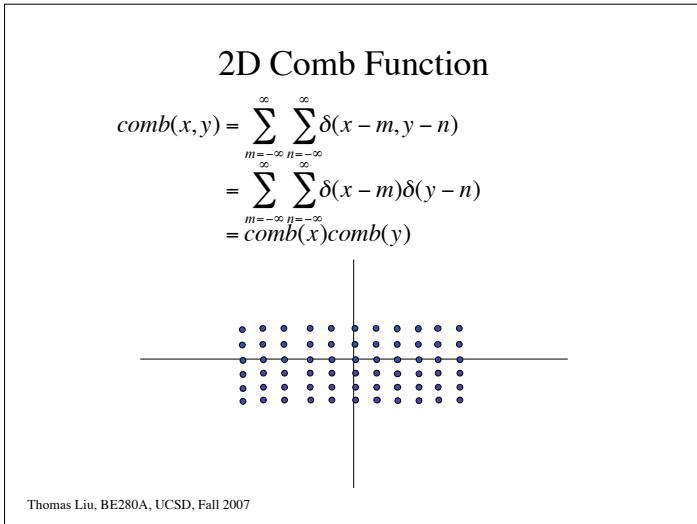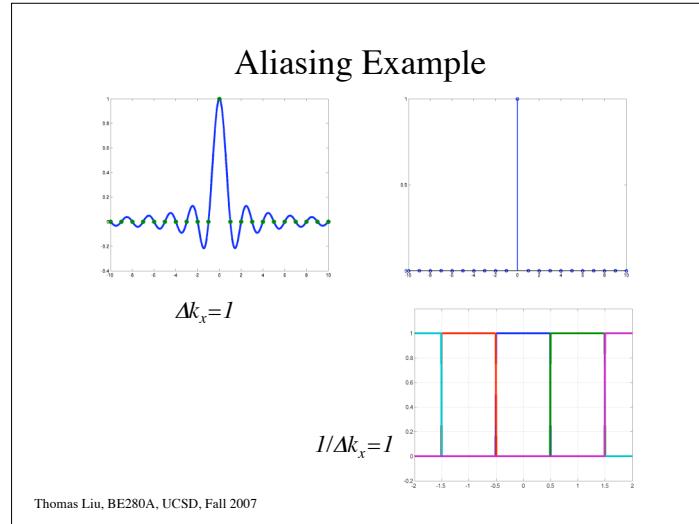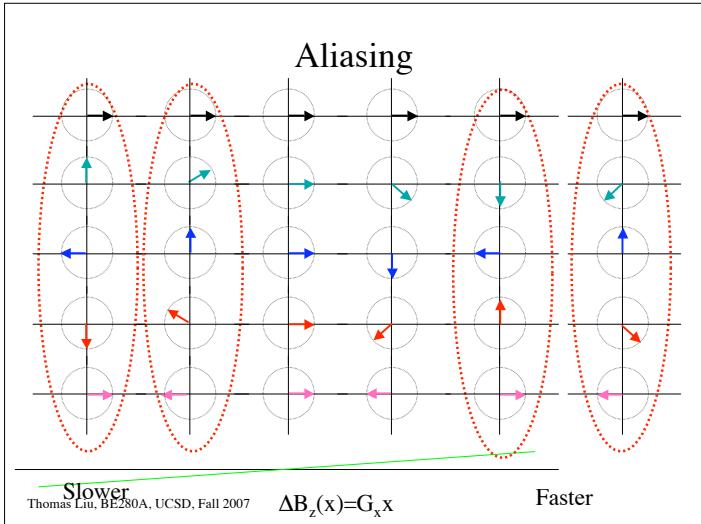

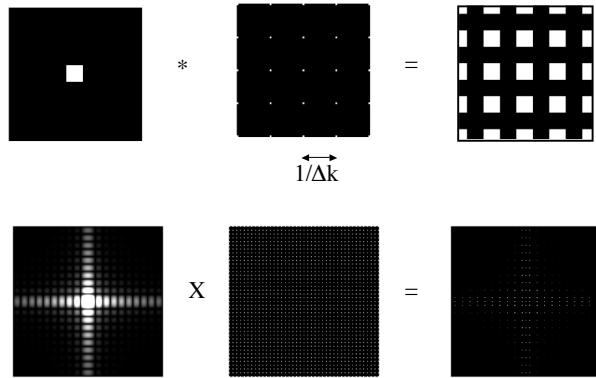

Thomas Liu, BE280A, UCSD, Fall 2007

## 2D k-space sampling

$$\begin{aligned}
 G_S(k_x, k_y) &= G(k_x, k_y) \frac{1}{\Delta k_x \Delta k_y} \text{comb}\left(\frac{k_x}{\Delta k_x}, \frac{k_y}{\Delta k_y}\right) \\
 &= G(k_x, k_y) \sum_{m=-\infty}^{\infty} \sum_{n=-\infty}^{\infty} \delta(k_x - m\Delta k_x, k_y - n\Delta k_y) \\
 &= \sum_{m=-\infty}^{\infty} \sum_{n=-\infty}^{\infty} G(m\Delta k_x, n\Delta k_y) \delta(k_x - m\Delta k_x, k_y - n\Delta k_y)
 \end{aligned}$$


Thomas Liu, BE280A, UCSD, Fall 2007

## 2D k-space sampling





$$\begin{aligned}
 g_s(x, y) &= F^{-1}[G_S(k_x, k_y)] \\
 &= F^{-1}\left[G(k_x, k_y) \frac{1}{\Delta k_x \Delta k_y} \text{comb}\left(\frac{k_x}{\Delta k_x}, \frac{k_y}{\Delta k_y}\right)\right] \\
 &= F^{-1}[G(k_x, k_y)] * F^{-1}\left[\frac{1}{\Delta k_x \Delta k_y} \text{comb}\left(\frac{k_x}{\Delta k_x}, \frac{k_y}{\Delta k_y}\right)\right] \\
 &= g(x, y) * * \text{comb}(x\Delta k_x) \text{comb}(y\Delta k_y) \\
 &= g(x) * * \sum_{m=-\infty}^{\infty} \sum_{n=-\infty}^{\infty} \delta(x\Delta k_x - m) \delta(y\Delta k_y - n) \\
 &= g(x) * * \frac{1}{\Delta k_x \Delta k_y} \sum_{m=-\infty}^{\infty} \sum_{n=-\infty}^{\infty} \delta(x - \frac{m}{\Delta k_x}) \delta(y - \frac{n}{\Delta k_y}) \\
 &= \frac{1}{\Delta k_x \Delta k_y} \sum_{m=-\infty}^{\infty} \sum_{n=-\infty}^{\infty} g(x - \frac{m}{\Delta k_x}, y - \frac{n}{\Delta k_y})
 \end{aligned}$$

Thomas Liu, BE280A, UCSD, Fall 2007

## Nyquist Conditions

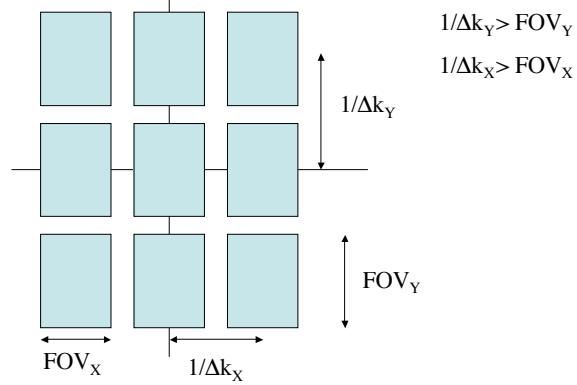


Thomas Liu, BE280A, UCSD, Fall 2007

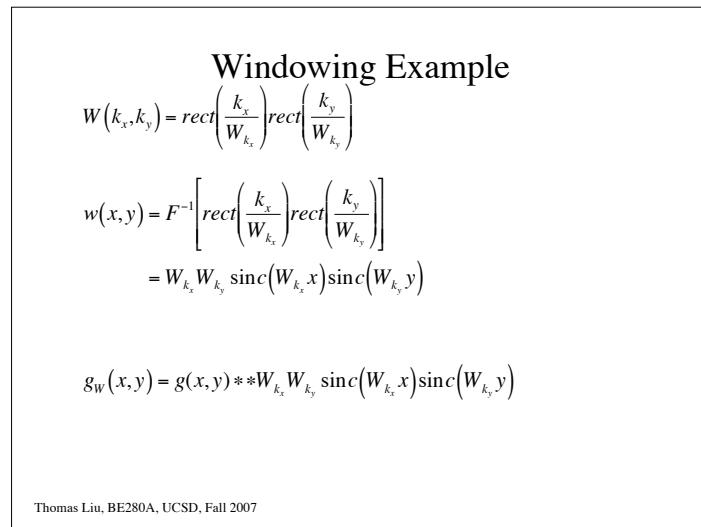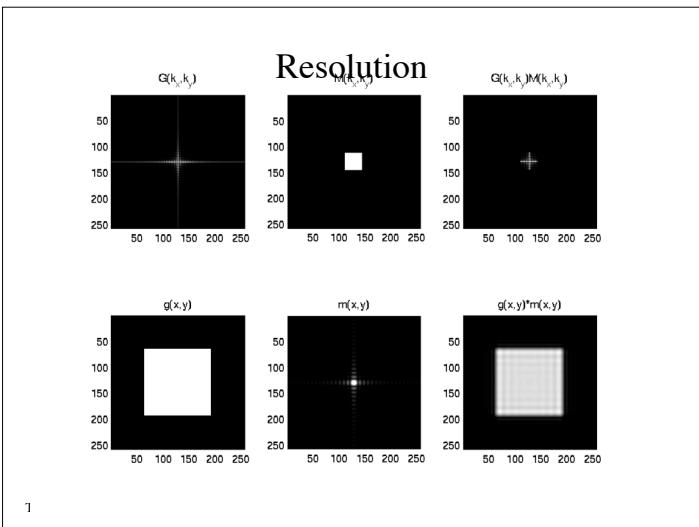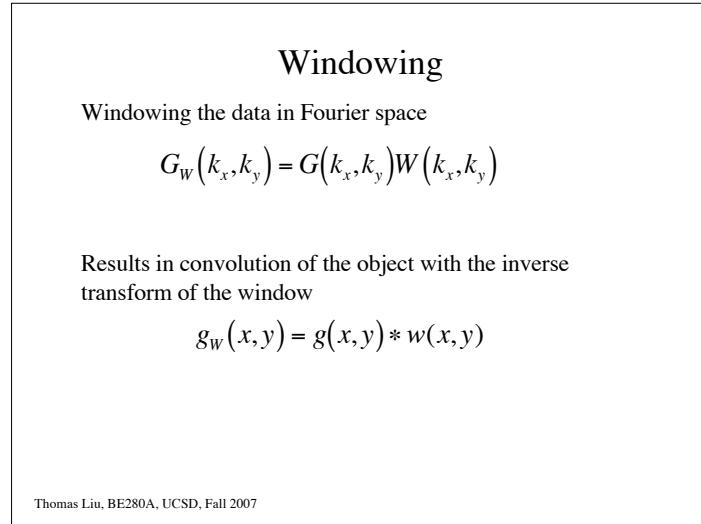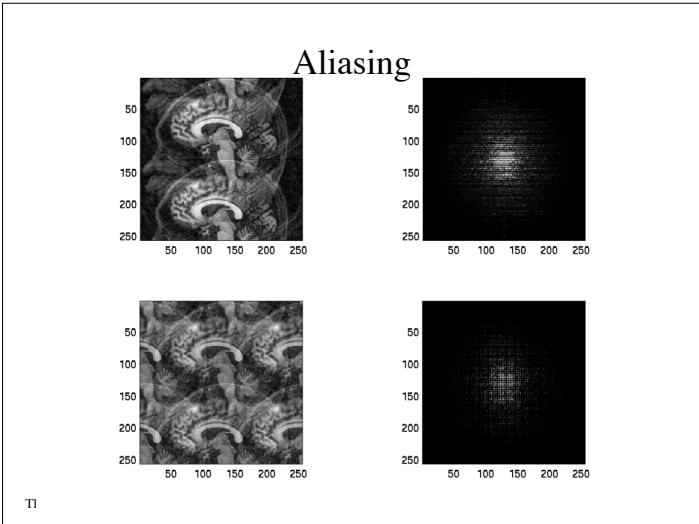


## Windowing Example

$$g(x, y) = [\delta(x) + \delta(x-1)]\delta(y)$$


$$\begin{aligned} g_w(x, y) &= [\delta(x) + \delta(x-1)]\delta(y) * W_{k_x} W_{k_y} \text{sinc}(W_{k_x}x)\text{sinc}(W_{k_y}y) \\ &= W_{k_x} W_{k_y} ([\delta(x) + \delta(x-1)] * \text{sinc}(W_{k_x}x))\text{sinc}(W_{k_y}y) \\ &= W_{k_x} W_{k_y} (\text{sinc}(W_{k_x}x) + \text{sinc}(W_{k_x}(x-1)))\text{sinc}(W_{k_y}y) \end{aligned}$$

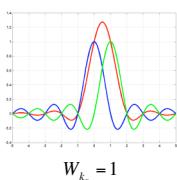

$$W_{k_x} = 1$$

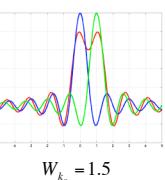

$$W_{k_x} = 1.5$$

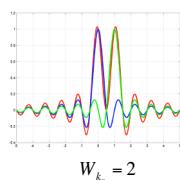

$$W_{k_x} = 2$$


Thomas Liu, BE280A, UCSD, Fall 2007

## Effective Width

$$w_E = \frac{1}{w(0)} \int_{-\infty}^{\infty} w(x) dx$$

Example

$$\begin{aligned} w_E &= \frac{1}{1} \int_{-\infty}^{\infty} \text{sinc}(W_{k_x}x) dx \\ &= F[\text{sinc}(W_{k_x}x)]|_{k_x=0} \\ &= \frac{1}{W_{k_x}} \text{rect}\left(\frac{k_x}{W_{k_x}}\right)|_{k_x=0} \\ &= \frac{1}{W_{k_x}} \end{aligned}$$

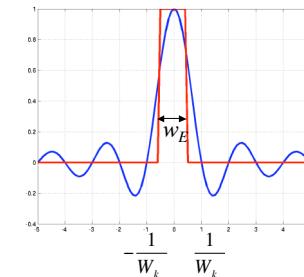



Thomas Liu, BE280A, UCSD, Fall 2007

## Resolution and spatial frequency

With a window of width  $W_{k_x}$ , the highest spatial frequency is  $W_{k_x}/2$ .

This corresponds to a spatial period of  $2/W_{k_x}$ .

$$\frac{1}{W_{k_x}} = \text{Effective Width}$$

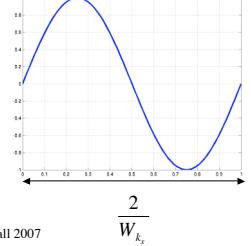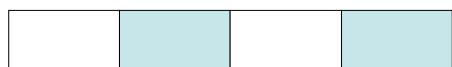

Thomas Liu, BE280A, UCSD, Fall 2007

## Sampling and Windowing

$$\begin{matrix} \text{Sampling} \\ \text{Matrix} \end{matrix} * \begin{matrix} \text{Window} \\ \text{Matrix} \end{matrix} * \begin{matrix} \text{Filter} \\ \text{Matrix} \end{matrix} = \begin{matrix} \text{Result} \\ \text{Matrix} \end{matrix}$$

$$\begin{matrix} \text{Sampling} \\ \text{Matrix} \end{matrix} X \begin{matrix} \text{Window} \\ \text{Matrix} \end{matrix} X \begin{matrix} \text{Filter} \\ \text{Matrix} \end{matrix} = \begin{matrix} \text{Result} \\ \text{Matrix} \end{matrix}$$

Thomas Liu, BE280A, UCSD, Fall 2007

## Sampling and Windowing

Sampling and windowing the data in Fourier space


$$G_{SW}(k_x, k_y) = G(k_x, k_y) \frac{1}{\Delta k_x \Delta k_y} \text{comb}\left(\frac{k_x}{\Delta k_x}, \frac{k_y}{\Delta k_y}\right) \text{rect}\left(\frac{k_x}{W_{k_x}}, \frac{k_y}{W_{k_y}}\right)$$

Results in replication and convolution in object space.

$$g_{SW}(x, y) = W_{k_x} W_{k_y} g(x, y) * * \text{comb}(\Delta k_x x, \Delta k_y y) * * \text{sinc}(W_{k_x} x) \text{sinc}(W_{k_y} y)$$


Thomas Liu, BE280A, UCSD, Fall 2007

## Sampling in $k_y$

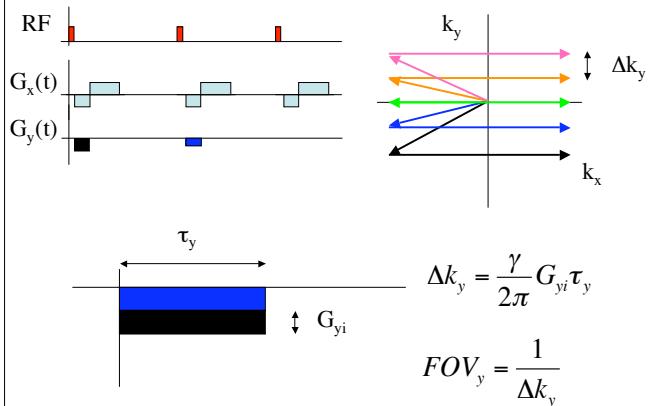

Thomas Liu, BE280A, UCSD, Fall 2007


## Sampling in $k_x$

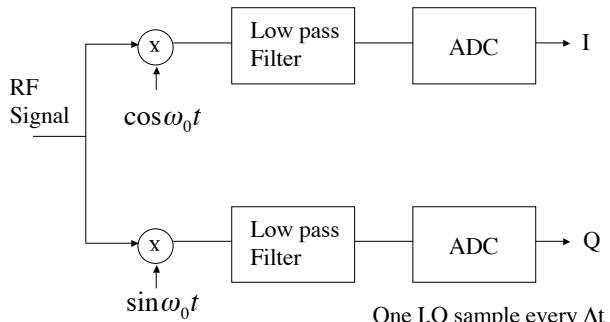





Note: In practice, there are number of ways of implementing this processing.

Thomas Liu, BE280A, UCSD, Fall 2007

## Sampling in $k_x$

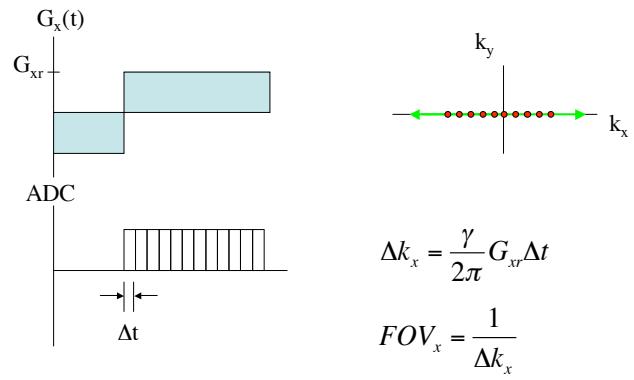

Thomas Liu, BE280A, UCSD, Fall 2007

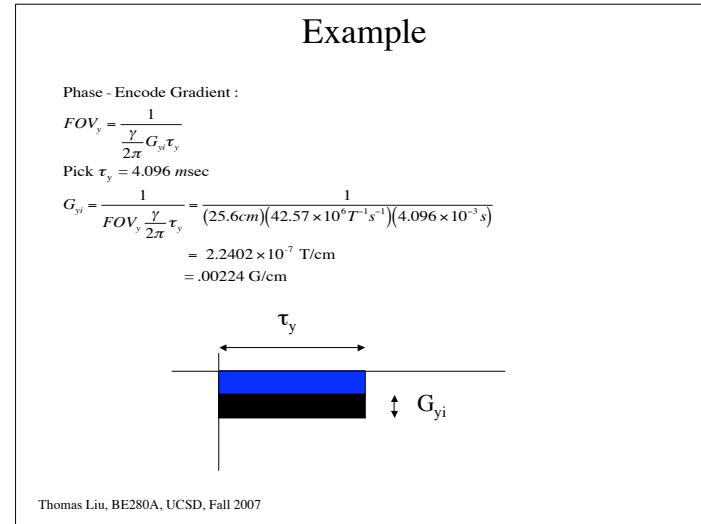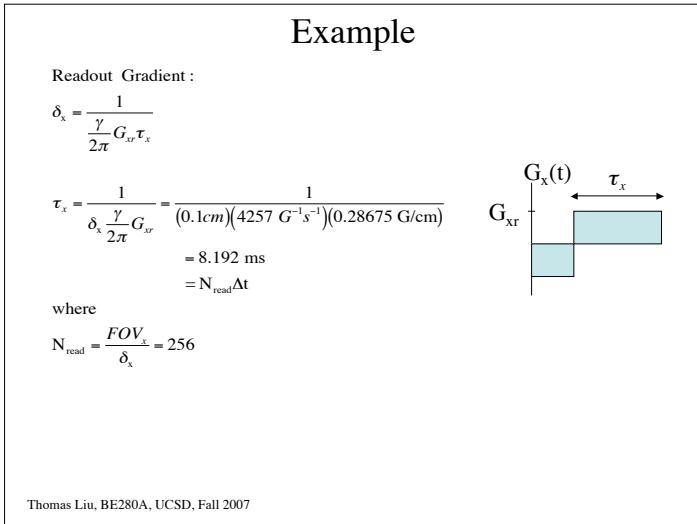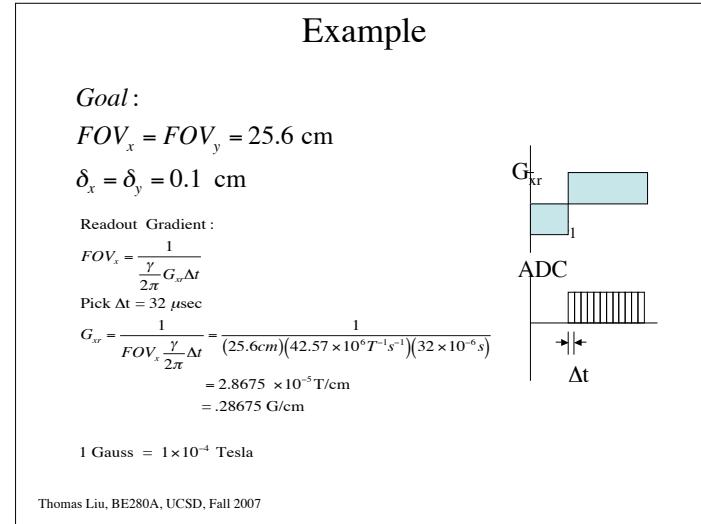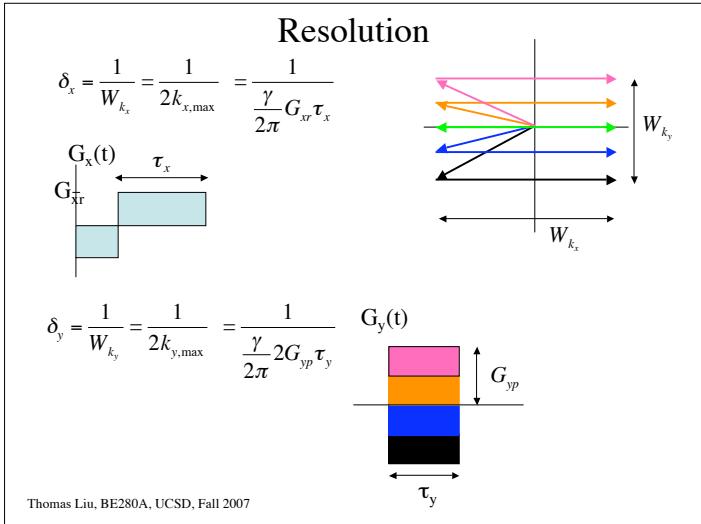

## Example


Phase - Encode Gradient :

$$\delta_y = \frac{1}{2\pi} \frac{\gamma}{2G_{yp}\tau_y}$$

$$G_{yp} = \frac{1}{\delta_y 2 \frac{\gamma}{2\pi} \tau_y} = \frac{1}{(0.1cm)(4257 G^{-1}s^{-1})(4.096 \times 10^{-3} s)} \\ = 0.2868 G/cm \\ = \frac{N_p}{2} G_{yi}$$

where


$$N_p = \frac{FOV_y}{\delta_y} = 256$$

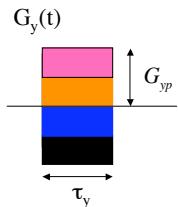

Thomas Liu, BE280A, UCSD, Fall 2007

## Sampling


In practice, an even number (typically power of 2) sample is usually taken in each direction to take advantage of the Fast Fourier Transform (FFT) for reconstruction.

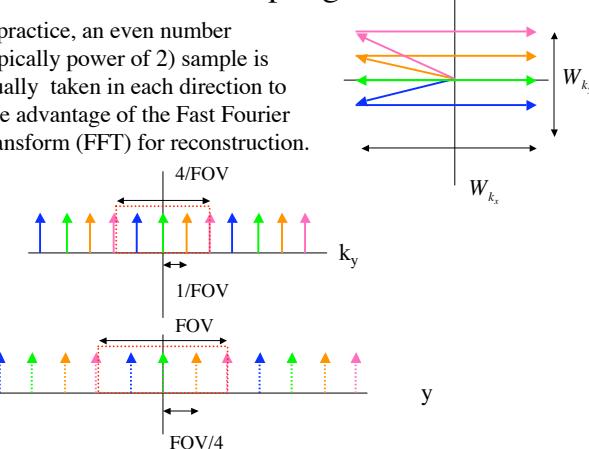

Thomas Liu, BE280A, UCSD, Fall 2007

## Example

Consider the k-space trajectory shown below. ADC samples are acquired at the points shown with  $\Delta t = 10 \mu\text{sec}$ . The desired FOV (both x and y) is 10 cm and the desired resolution (both x and y) is 2.5 cm. Draw the gradient waveforms required to achieve the k-space trajectory. Label the waveform with the gradient amplitudes required to achieve the desired FOV and resolution. Also, make sure to label the time axis correctly.

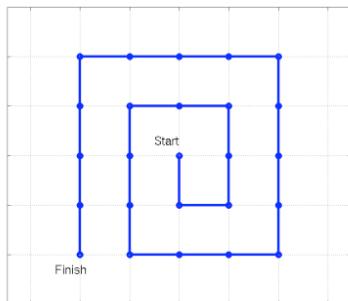

Thomas Liu, BE280A, UCSD, Fall 2007

## SCAN TIMING

# of Echoes	1	2	4
TE Min Full	1	2	4
TE	1	2	4
TR	750	1500	3000
Inv. Time	1	2	4
T12	1	2	4
Flip Angle	1	2	4
Echo Train Length	1	2	4
Bandwidth	2.5	5.0	10.0
Bandwidth2	1	2	4

## ACQUISITION TIMING





Freq	352	Freq DIR	A/P
Phase	192	Center Freq	Water
NEX	2.0	Flow Comp	Direction
Phase FOV	0.75	Autoshim	Phase Correct
# of Acqs Before Pause	1	Contrast Amt	0.0
Agent	None		

## SCANNING RANGE


FOV	22	S/I	L/R Center	P/A Center
Slice	5.0	Start	End	# Slices
Spacing	2.0	Table Delta	Actual End	Table Delta

Thomas Liu, BE280A, UCSD, Fall 2007


GE Medical Systems 2003

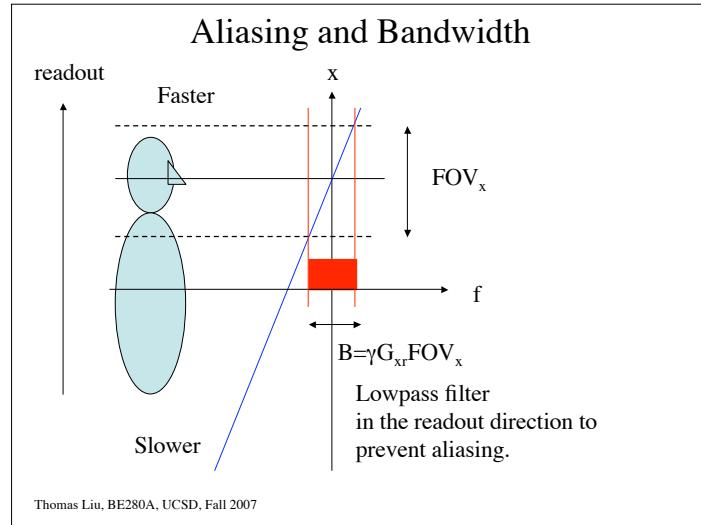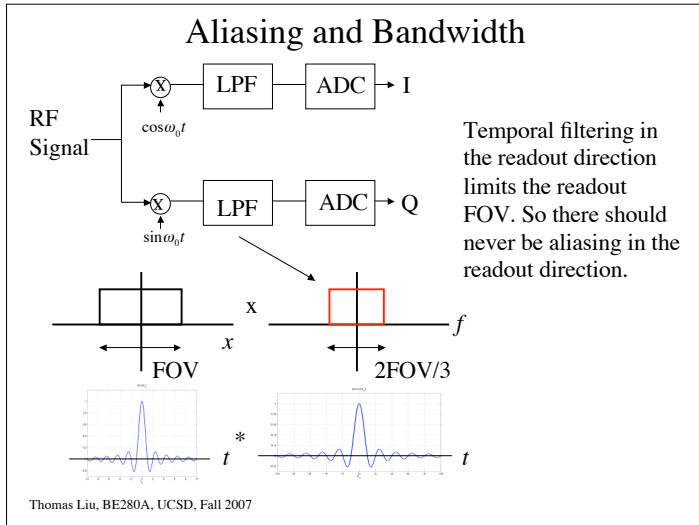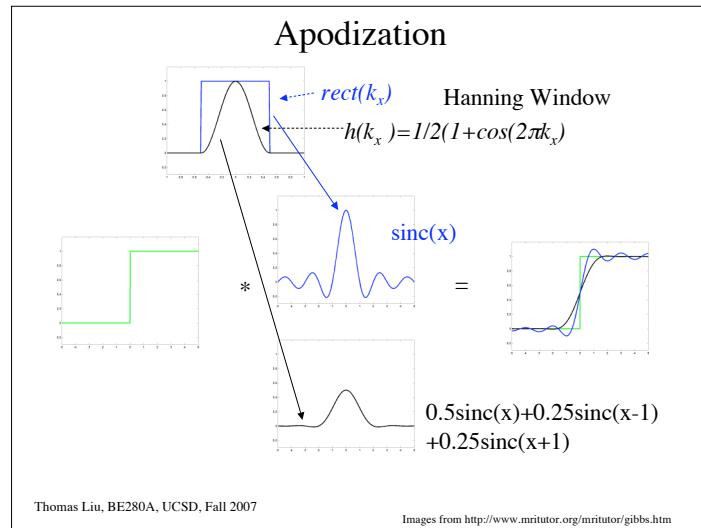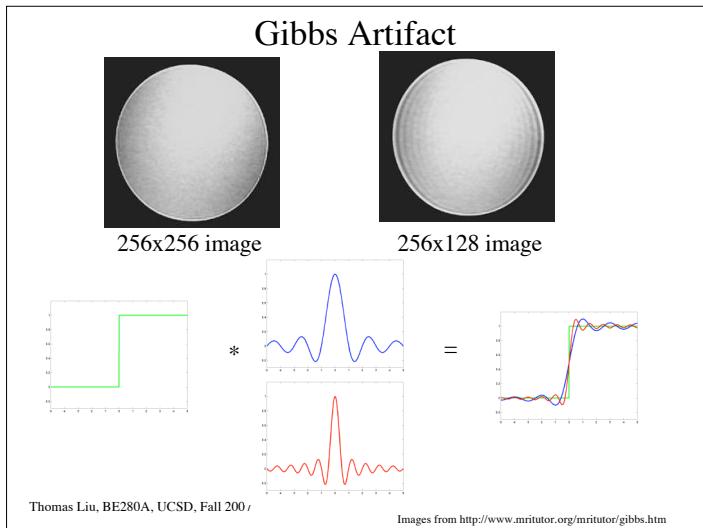

**Figure 7-31 Default Axial Directions**

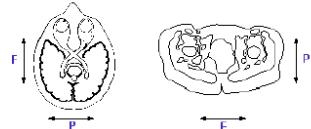

**Figure 7-32 Default Sagittal and Coronal Directions**

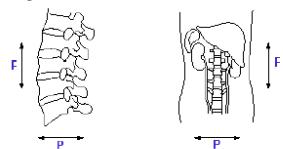

Thomas Liu, BE280A, UCSD, Fall 2007

GE Medical Systems 2003